## Digitale Personenwaage Linette/Lilian

Digital bathroom scale Linette/Lilian | Báscula de baño digital Linette/Lilian | Pèse-personne numérique Linette/Lilian | Bilancia digitale pesapersona Linette/Lilian | Cyfrowa waga osobowa Linette/Lilian

Deutsch

English

Español

Français

Italiano

Polski

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BE 1906 | 1907

Operating instructions | Instrucciones de operación | Mode d'emploi  
Istruzioni per l'uso | Instrukcja obsługi

# Liebe Kundin, lieber Kunde !

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# Dear Customer!

You have chosen to purchase a high-quality ADE brand product , which combines intelligent functions with an exceptional design. These bathroom scales make it possible for you to determine and monitor your body weight. The wealth of experience of the ADE brand ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

## General information

---

### About this operating manual

This operating manual describes how to operate and maintain your product safely. Keep this operating manual in a safe place for future reference. If you pass this product on to someone else, pass on this operating manual as well.

Failure to observe this operating manual may result in injuries or damage to your product.

### Explanation of symbols

This symbol combined with the word DANGER warns of serious injuries.

This symbol combined with the word WARNING warns of moderate and mild injuries.

This symbol combined with the word NOTE warns of material damage.

This symbol indicates additional information and general notes.

# Table of Contents

---

General information .....	12
Technical data.....	13
Intended use .....	14
Safety .....	14
Overview .....	16
Scope of delivery .....	16
Starting up .....	17
Weighing .....	18
Body mass index (BMI) .....	18
Resetting the scale.....	19
Cleaning.....	20
Fault/Remedy.....	20
Declaration of Conformity.....	21
Warranty .....	21
Disposal.....	21

## Technical data

---

Model:	Digital bathroom scale BE 1906 "Linette" (black), BE 1907 "Lilian" (silver)
Battery:	2x LR03 (AAA)/1.5 V 
Rated current:	BE 1906: 15 mA BE 1907: 10 mA
Measuring range:	Max. 180 kg
Variance:	1 % +/- 100g
Display scaling:	50 g
Scale dimensions:	Approx. 333 x 303 x 23 mm
Distributor (no service address):	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

## Intended use

---

The bathroom scale is intended for private use, for weighing persons and determining body mass index. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

## Safety

---

- This product may be used by children aged 8 and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge if they are supervised or have been instructed in the safe use of the product and understand the risks that may arise.
- Children must not play with the product.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. You should therefore keep the batteries and scale out of reach of small children. If a battery has been swallowed, medical assistance must be sought immediately.

## DANGER for children

- Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

## RISK of fire/burning and/or explosion

- There is a risk of explosion if the batteries are replaced improperly. Replace the battery only with the same or a similar type. Pay attention to the "Technical data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise, there is an increased risk of explosion and leakage; furthermore, gases may escape!
- Avoid contact of the battery fluid with skin, eyes, and mucosa. In case of contact, immediately rinse the affected areas with plenty of clean water and consult a doctor.

## **WARNING of risk of injury**


- Do not step onto the outer edge of any one side of the scales. Risk of tipping!
- Do not step onto the scales with wet feet or when its surface is damp. Risk of slipping!
- If the scale falls and shatters, there is a risk of injury from glass splinters.

## **BEWARE of material damage**

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the “Technical data”. Place yourself on the scales exclusively for weighing purposes.
- Protect the scale from strong shocks and vibrations.
- Remove the battery from the scale when it is depleted or when you are no longer using the scale. You will thus avoid damage that may result from leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not subject the battery to any extreme conditions, by – for example – storing it on radiators or under direct sunlight.  
Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the product. Only allow repairs to be carried out by a specialist workshop. Repairs carried out improperly can considerably endanger the user.
- Do not place the scales in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indicators or malfunctions.
- There are feet on the underside of the scale. Since surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place a non-slip pad under the scale.

# Overview

---

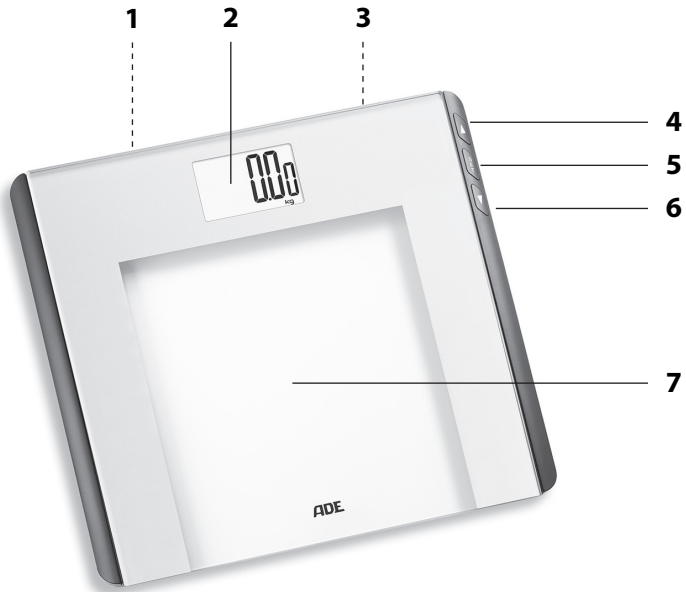- 1 UNIT button (underside): Changing the unit of measurement
- 2 Display
- 3 Battery compartment (underside)
- 4 Button ▲: Select memory cell, increase data
- 5 SET button: Bring up input mode, save input
- 6 Button ▼: Select memory cell, reduce data
- 7 Weighing surface

## Scope of delivery

---

- Bathroom Scales
- Battery type LR03 (AAA), 1.5 V  (2x)
- Operating Manual

# Starting up

## Inserting/changing batteries

For use, you need 2 batteries of type LR03 (AAA), 1.5 V (included). The battery compartment is located on the underside of the device.

1. Push the locking tab towards the edge of the scale and remove the battery compartment cover.
2. Insert the batteries as indicated at the bottom of the battery compartment. Pay attention to the correct polarity (+/-).
3. Replace the battery compartment cover. The locking tab must snap audibly into place.

Suggestions for changing the battery:

- A battery change is necessary when "Lo" or nothing at all is shown on the display.
- Only use the type of battery that is specified in the "Technical data".
- Always change all batteries at the same time.
- Dispose of old batteries in an environmentally appropriate way (see "Disposal").
- Calibrate the scale, see "Setting the scale".

## Choosing the unit of measurement

When the scales are delivered, the unit of measurement is set to "kg" as standard. However, the measurement result can also be displayed in the English unit "lb" (pound) or "st lb" (stone:pound).

1. To switch on the scale, step on the middle of the weighing surface.
2. Press repeatedly the Unit button on the back of the scale.  
Note the display.

## Positioning the scales

1. Place the scale on an even, firm and dry surface. Do not place the scales on a carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows "0.0", calibration is complete.  
This calibration must be performed after each battery change.

# Weighing

---

## **WARNING** against injuries

- Do not step onto the scales with wet feet or when its surface is damp. Risk of slipping!
- Position yourself on the scales. Remain standing still and do not hold on to anything.  
The measured weight will be displayed shortly thereafter.  
A few seconds after you step off the weighing surface, the scales will switch off automatically.

- Weigh yourself with as little clothing on as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

## **Body mass index (BMI)**

---

The scale can determine your personal body mass index (BMI). For more information, see section "What your body mass index means".

### **Entering your height**

To determine your BMI, you must enter your height in centimetres once and save it. The scale offers a total of 12 memory cells.

1. Press and hold SET for a few seconds.  
A memory cell number ("P-01") will be shown blinking on the display.
2. Using the ▲ or the ▼ button, select the desired memory cell.
3. Press and hold SET for a few seconds.  
The pre-set value for your height in centimetres will blink on the display.
4. Enter your height by repeatedly pressing the ▲ or ▼ button.
5. Finally, press SET.

No further input is required to determine your BMI.

If you do nothing else, the scales will automatically switch off after a few seconds.

## Determining weight and BMI

1. Press and hold SET for a few seconds.  
A memory cell number will be shown blinking on the display.
2. Select the memory cell in which your height is saved using the ▲ or ▼ button. Wait till "0.0" is shown on the display.
3. Position yourself on the scales. Remain standing still and do not hold on to anything.
  - First the measured weight is displayed.
  - Then the memory cell is displayed along with the measured weight.
  - Subsequently, you will see your body mass index displayed (see the next section).

The weighing result (weight and BMI) is displayed three times in total. The scale then automatically switches off.

## What your body mass index means

The BMI indicates the ratio of height to weight and is only an approximate standard value. According to the World Health Organization (WHO), an index of 18.6 to 25 is optimal.

Since men typically have a higher percentage of muscle mass in the total body mass than women do, the age and gender play an important role in interpreting the BMI.

BMI (kg/m <sup>2</sup> )	Body weight
< 18.5	Underweight
18.6 - 25	Normal weight
25.1 - 30	Overweight
30.1 and above	Obesity

## Resetting the scale

Reset the scale to its default settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment and wait a moment.
2. Reinsert the batteries.  
Pay attention to the correct polarity of the batteries (+/-).

After switching the scales on, you can weigh yourself as usual.

# Cleaning

## NOTE on material damage

- Do not immerse the scale in water or other liquids.
  - Do not under any circumstances use harsh, scouring or abrasive cleaning agents or the like to clean the scale. Otherwise, the surface could be scratched.
- 
- Wipe the scale as needed with a slightly damp soft cloth.
  - Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

# Fault/Remedy

Fault	Possible cause
No function	<ul style="list-style-type: none"><li>▪ Are the batteries empty or inserted the wrong way around?</li></ul>
The scale appears to be displaying wrong weighing results	<ul style="list-style-type: none"><li>▪ Is the scale located on a firm, even surface?</li><li>▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale?</li><li>▪ Remove the batteries from the scale and wait a moment. Reinsert the battery.</li></ul>
“Err” is shown on the display	<ul style="list-style-type: none"><li>▪ The maximum capacity of the scale (180 kg) has been exceeded.</li></ul>
“Lo” is shown on the display	<ul style="list-style-type: none"><li>▪ The batteries are empty (see “Inserting/ changing batteries”).</li></ul>
“C” is shown on the display	<ul style="list-style-type: none"><li>▪ The weighing/measurement process was faulty: The scales are not standing stably or you are not standing still on the weighing surface during weighing. Repeat the process.</li></ul>

# Declaration of Conformity

Waagen-Schmitt GmbH hereby declares that the electronic bathroom scale BE 1906/BE 1907 complies with EU directives 2014/30/EU and 2011/65/EU. This declaration will lose its validity if any modification is made to the scale without our consent.

Hamburg, March 2019  
Waagen-Schmitt GmbH

## Warranty

Waagen-Schmitt GmbH warrants for 2 years from the purchase date the cost-free rectification of defects in material and workmanship through repair or exchange. In the event of a warranty claim, please return the product with the proof of purchase (specifying the reason for claim) to your dealer.

## Disposal

### Disposing of the packaging

Dispose of the packaging with similar materials. Dispose of cardboard with waste paper and of films with recyclable materials.

### Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

#### **Devices must not be disposed of in normal household waste.**

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Bring the old device to a collection point for electronic waste or to a recycling depot. For more information, contact your local waste disposal firm or your local administrative body.

#### **Standard and rechargeable batteries should not be disposed of along with household waste.**

All standard and rechargeable batteries must be brought to a local or district collection point or to the retailer. Such standard and rechargeable batteries can thus be disposed of in an environmentally friendly way.

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